




SILVERWAREPOS


HOW TO CALIBRATE HP TOUCHSCREEN (Windows 7)

This document will walk you through how to calibrate an HP workstation using Windows 7.

METHOD 1: Using the Microsoft calibration tool


To calibrate the screen of the HP TouchSmart using the Microsoft calibration tool, complete the following steps.

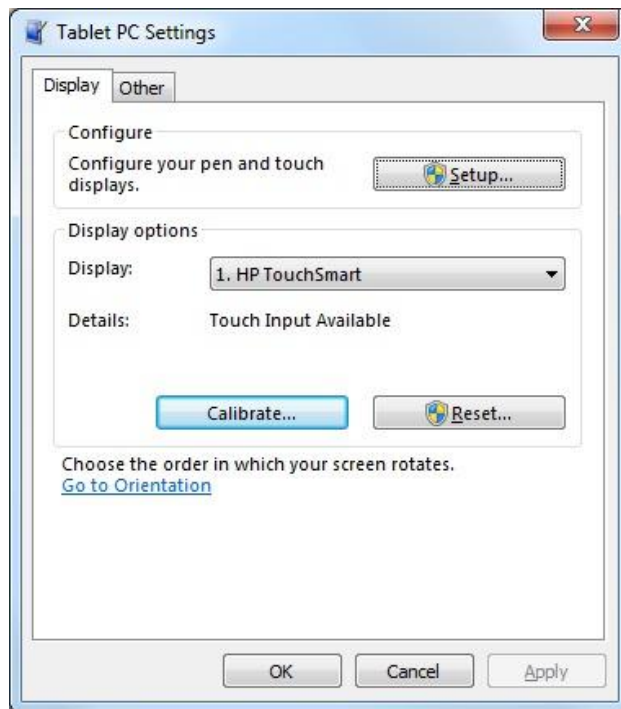1. Clean the screen before calibrating.
 - a. Turn the computer OFF and use a soft, damp cloth and gently wipe the screen to remove dirt and smudges. Give special attention to the upper left and right corners of the screen. Turn on the computer when finished.
2. Log into Windows Support Administrator.
3. Click or tap Start  on the task bar, and then click or tap Control Panel.
4. In Control Panel, tap Hardware and Sound, and then tap Tablet PC Settings.
5. On the Display tab, tap Calibrate.




6. Tap the crosshair with your finger each time that it displays on the screen.

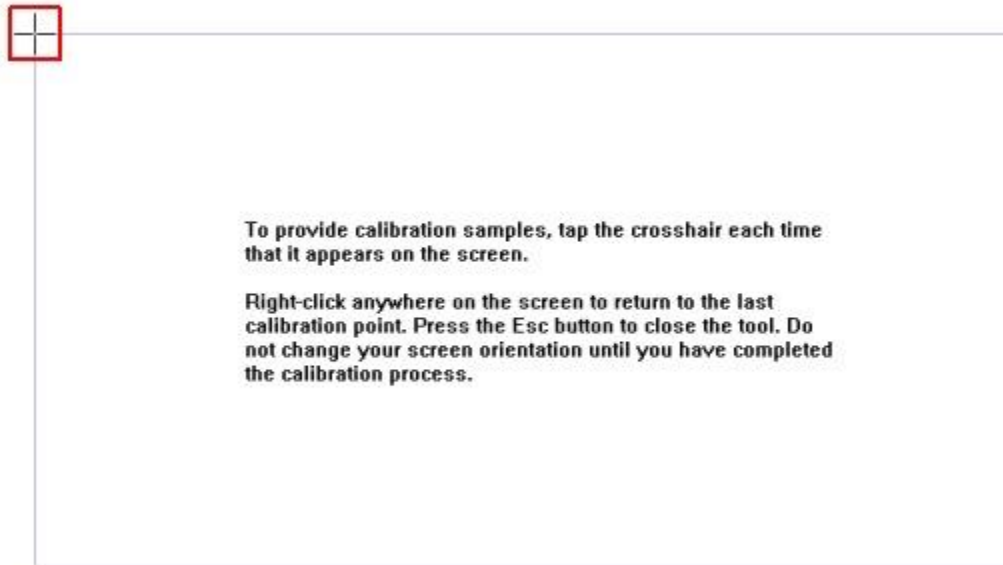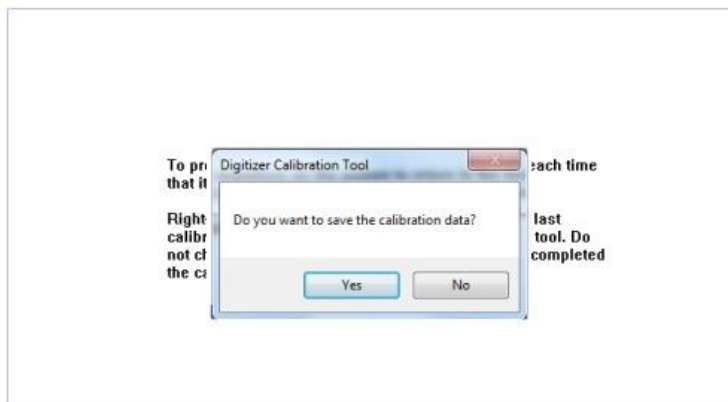NOTE: Tap the cross hairs quickly. Do not hold and press your finger or the stylus on the screen for a long time, or the tool might not display the next calibration point.




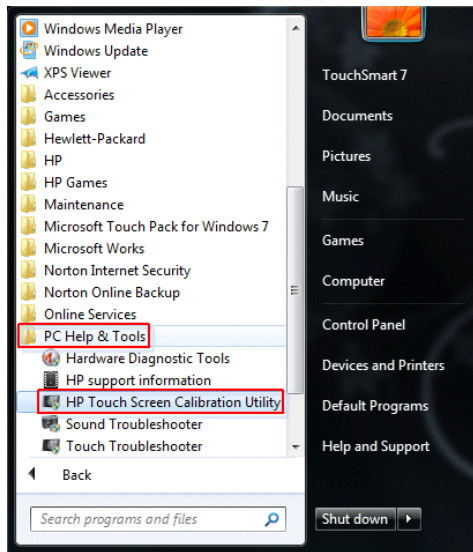
7. When done, tap Yes to save the calibration data.



8. If the screen is still not calibrated correctly, try resetting the computer back to the original calibration settings. Check the Reset button on the Display tab in the Tablet PC Settings.

METHOD 2: Using the HP calibration tool

1. Clean the screen before calibrating.
 - a. Turn the computer OFF and use a soft, damp cloth and gently wipe the screen to remove dirt and smudges. Give special attention to the upper left and right corners of the screen. Turn on the computer when finished.
2. Log into Windows Support Administrator.
3. Click or tap Start  on the task bar, and then click or tap All Programs.
4. Click or tap PC Help & Tools, and then click or tap HP Touch Screen Calibration Utility.



5. In the Calibration Utility, the cursor displays in the upper left corner of the screen. Touch the cursor for two seconds and release. The cursor moves to the next corner. Continue touching and releasing the cursor in each corner.
6. After the fourth corner, click the Close button.

